

Accept No Substitutes

Durcon's time-tested Epoxy Resin worktop is the industry standard for harsh laboratory environments. It is safe, durable and widely available for projects of all types and sizes. Durcon offers several edge and backsplash styles to meet various requirements for safety and comfort. Each style can be configured into virtually any furniture configuration and can be installed on permanent or freestanding casework, tables and carts.

ClassicTop Bevel Edge




One inch thick Black Onyx ClassicTops with machine bevelled edges are the laboratory industry standard. Available in seamless lengths up to 96" [2438mm], ClassicTops provide the greatest flexibility for use in every application.

Upon request, 5/8", 3/4" or 1-1/4" [15, 19 or 32mm] thick surfaces in standard and custom colors are available.

ClassicTop Radius Edge

ClassicTops with a 1/4" machined radius feature the latest worksurface edge option from Durcon. These worksurfaces are molded and fabricated to the same exacting standards as our machine bevelled ClassicTops, but feature machine radiused edges and eased corners which can provide additional safety and comfort for laboratory users.



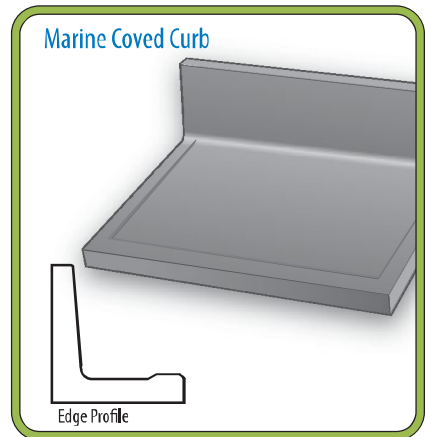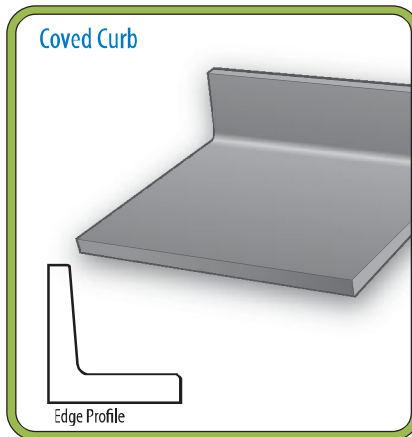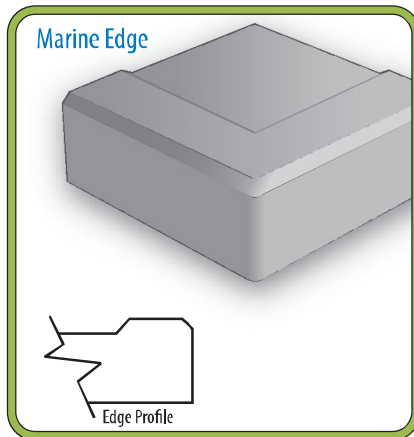

Classic MarineTop

Classic MarineTop worksurfaces are recommended for areas where liquid containment is a priority. These 1/4" [6mm] dished worksurfaces help protect casework finishes, flooring and lab personnel by keeping most spills on the countertop.

ClassicTop Coved Curb

Classic Coved Curb worksurfaces should be used where eliminating backsplash joints is the highest priority. They are available in seamless lengths up to 98" [2489mm] and depths up to 38" [965mm]. Classic Coved Curb worksurfaces are molded 1" [25mm] thick. By special request, 3/4" or 1-1/4" [19mm or 32mm] thick worksurfaces can be ordered. Integral curbs are 4" [102mm] tall.

For special curb height requirements, contact our customer service department.

Classic Marine Coved Curb

Classic Marine Coved Curb worksurfaces are recommended for use where liquid containment and eliminating backsplash joints is the highest priority. The 1/4" [6mm] deep dish helps keep liquids from overflowing. Marine Coved Curb worksurfaces are either 1" or 1-1/4" [25 or 32mm] thick at the rail. The standard 4" [102mm] high integral curbs maintain a uniform backsplash height when installed in-line with non-marine edged Classic Coved Curb worksurfaces.

