


LabArt™ Pegboard

Features:

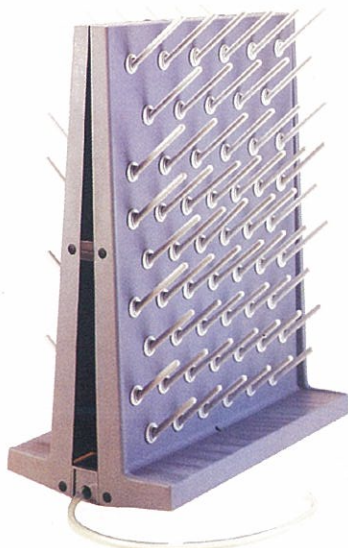

1. Injection molded main structure provides high rigidity and lighter weight.
2. Recycled polypropylene (PP) pegboard is easy to clean and environment friendly.
3. Pegs fit snugly into molded insertion holes that are specially designed to keep pegs in a uniform upright 40° angle for stable glassware storage.
4. Removable pegs, of unique tenon design, are plugged to the pegboard with convenient front-side mounting/demounting mechanism.
5. Fifty-five peg holes for each pegboard, special water stopping design to prevent spill.
6. Removable pegs are easily relocated to accommodate glassware items of various sizes, achieving high mobility and flexibility. Un-used holes are covered with caps for esthetics.
7. The tilted drain groove design facilitates easy and thorough drainage.
8. Two types of pegboard, single side for wall mounting and double sides for bench mounting, are available.
9. The loading capacity of each peg is up to 7 kg (15.4 lbs).
10. Honored USA and China patents.

Standard Accessories:

Description	Quantity	Unit	Required Quantity For Wall Mounting	Required Quantity For Bench Mounting
100 mm Peg	20	pc	19	38
120 mm Peg	40	pc	36	72
Plug (For empty holes only)	28	pc	For empty holes only	For empty holes only
Upper Bracket	2	pc	2	2
Middle Bracket	1	pc	0	2
Lower Mounting Kit	1	pc	0	2
Lower Bracket	1	pc	1	0
L Type Connector	1	pc	1	2
Y Type Connector	1	pc	0	1
Plastic Rivet	6	pc	4	12
Screw and Plastic Sleeve	4	set	4	2
Drain Hose, 1 meter	1	pc	1	1

Specification:

Dimension H x W x D	700x550x30/125mm
Color for pegboard	Platinum Gray, Black
Color for peg	Pearl Gray
Peg length	120, 100mm
Peg diameter	Φ 10mm
Material	Polypropylene (PP)
Insert holes for pegs	55 holes per side

